



# Digital Camera Adaptor TD-2

S e r v i n g   y o u r   v i s i o n



# 1 Light path changeable in accordance with right or left better eye through one-touch operation

The light path to capture images into the digital camera can be changed in accordance with the examiner's better eye. By adjusting the dial on the back, you can change the light path through one-touch operation. The fundus can be observed without worrying about displacement between an image in the field of view and a captured image.



# 2 Yellow filter IN/OUT mechanism


High-contrast images can be obtained by using a built-in yellow filter in cornea fluorescence observation. Easy filter IN or OUT is available through one-touch operation. In its setting with our slitlamp, high-contrast images can be obtained by the best matching with a cobalt blue filter.



# 3 Capturing of still image with one push

By connecting a trigger cable between the slitlamp and camera, a still image can be captured every time the joystick button of the slitlamp is pushed.

\*The trigger cable may not be available depending on the type of the slitlamp. For more information, contact us.

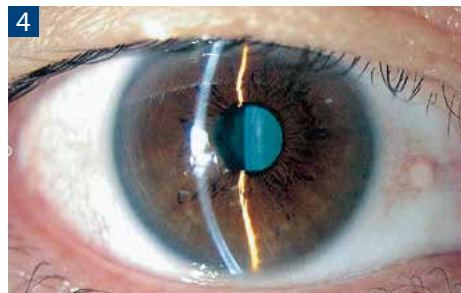

# 4 Combination with background illumination system

Combination with background illumination of the slitlamp is more effective in the explanation to patient.

\*Background illumination systems differ with slitlamps. For more information, contact us.



# 5 Slitlamp models compatible with TD-2 types

Compatible models	Type A	700GL, 4ZL, 2ZL, SM-70N, SM-90N, 300XL
-------------------	--------	--

\*For more information about 300XL, contact us.

# 6 Camera models compatible with TD-2

TD-2 types do not come with a digital camera. For the purchase of a digital camera, refer to the following compatible cameras:

Compatible models	CANON Digital SLR Cameras
	Nikon Digital SLR Cameras

\*For more information on compatible cameras, contact us.



● Design and specifications are subject to change as improvements are made to the product.